

**PAJUNK®**

## *TrocaSys*

*Extended field of application  
for instruments up to 12 mm*



Laparoscopy

TrocaSys **NEW**

# The disposable trocar system for instruments with a diameter of 5–12 mm

PAJUNK® expands its disposable trocar system with an additional access system with universal seal. This permits the gas-tight insertion of instruments with a diameter of 5 to 12 mm. The choice of instruments that can be used as well as the diameter increases accordingly.



**Lockable obturator**  
In completely inserted position, a bayonet lock connects the obturator securely to the valve housing.

**Simple unlocking**  
A twist is all that's needed.



**Universal seal**  
All instruments with diameters between 5 and 12 mm can be inserted completely gas-tight.

**Trocar and obturator**




TrocaSys PROTECT – atraumatic dilating trocar



Blunt obturator for the "Hasson technique"



Double seal system

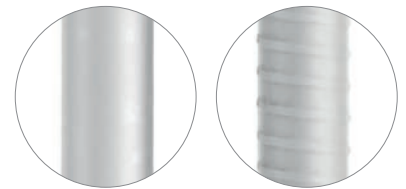

**Versions available:**  
**Smooth or corrugated guidance tube**  
The groove structure of the corrugated guide tube ensures that the trocar sleeve is securely held in the abdominal wall when the instruments are inserted and removed.

## TrocaSys



**TrocaSys PROTECT**  
Disposable trocar system with dilating trocar

Size	smooth guidance tube	corrugated guidance tube	PU
ø 12 mm	1289-03-12	1289-13-12	5



**TrocaSys**  
Disposable trocar system with blunt obturator

Size	smooth guidance tube	corrugated guidance tube	PU
ø 12 mm	1289-02-12	1289-12-12	5



**TrocaSys trocar sleeve single, without trocar / obturator**

Size	smooth	corrugated	PU
ø 12 mm	1289-09-12	1289-19-12	5

**PAJUNK Medical Systems L.P.**  
6611 Bay Circle, Suite 140  
Norcross, Georgia 30071  
Phone: 770.493.6832  
Fax: 678.514.3388  
E-mail: [info@pajunk-usa.com](mailto:info@pajunk-usa.com)